




# PATIENT KIT AND GUIDE

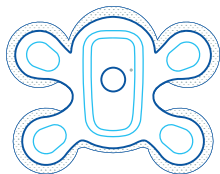

**2AYe**

Wearable Biosensor



# HOLTER MONITORING MADE EASY

You have been prescribed a Holter Test to continuously record your ECG, or your heart's electrical activity. Your heart data will automatically be transferred from the Biosensor, through your smartphone running the UbiqVue Holter App, to a secure server, where the data is analyzed, and a report created for your physician to review.

## WHAT'S COVERED:

A



### APPLY THE BIOSENSOR

Simple steps to apply the Biosensor

B



### CONNECT THE BIOSENSOR

Download the App and pair with Biosensor


C



### REMOVE THE BIOSENSOR

At the end of the procedure, remove your Biosensor

## WATCH INSTRUCTIONAL VIDEO

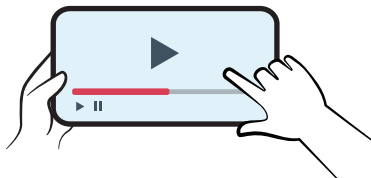

Scan this QR code to watch a video of these instructions

OR




Q 2AYe UbiqVue Patient Kit

Search for "2AYe UbiqVue Patient Kit" on YouTube



## BIOSENSOR APPLICATION

1

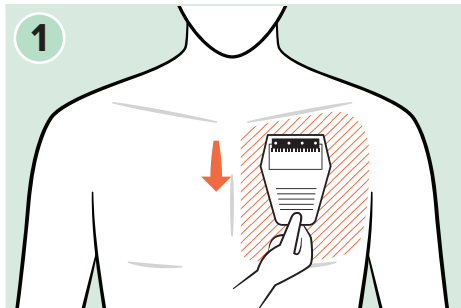


### SHAVE CHEST HAIR

If required, shave entire left side of your chest using the razor provided.



2

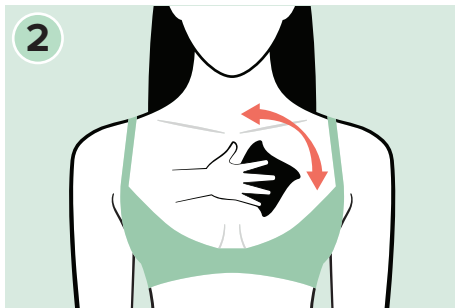


### PREPARE SKIN

Use the disinfectant wipe to clean the chest area.



3

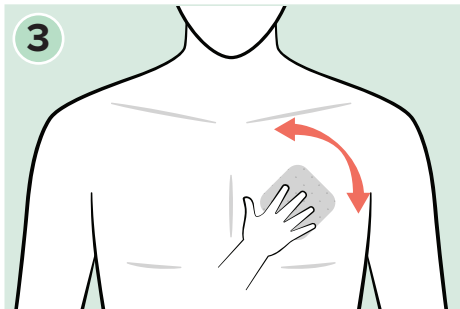


### PREPARE SKIN

Dry the chest area with the cotton pad




4

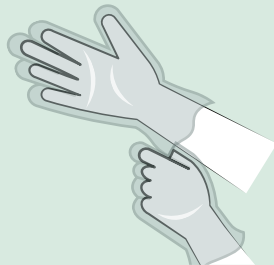

### APPLY ADHESIVE


Put on the gloves provided before using the sticky adhesive wipe



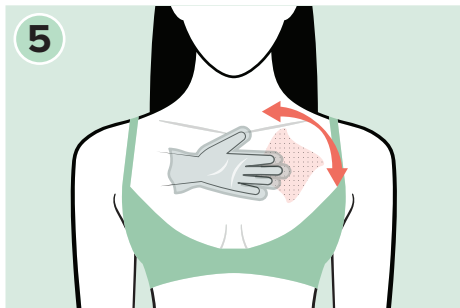

## APPLY ADHESIVE


Rub the adhesive wipe on the chest area and wait until the skin is completely dry

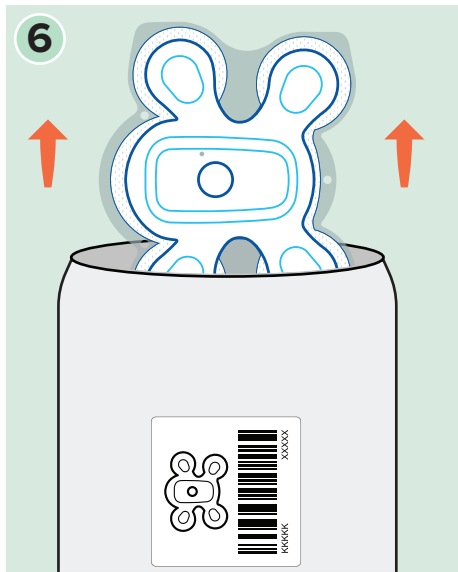

## UNPACK BIOSENSOR

Remove the Biosensor from the silver pouch



7

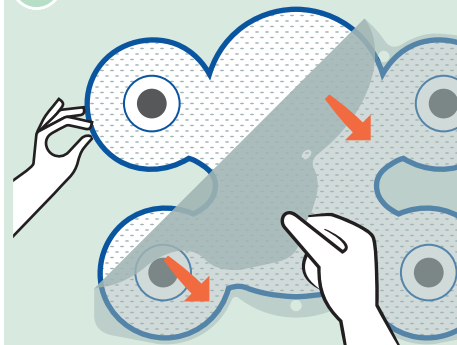

### PEEL BACKING FILM

Slowly peel off the Biosensor backing film, starting from one corner.



8

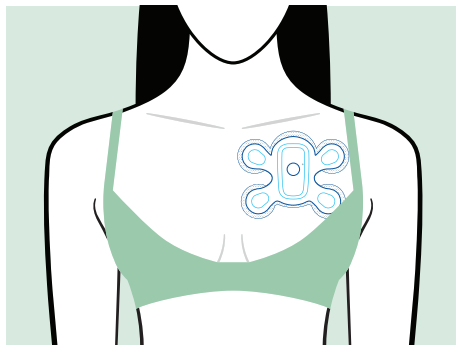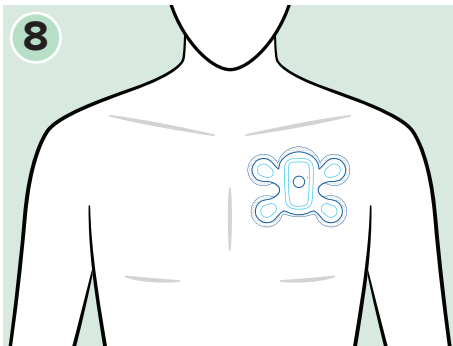

### APPLY BIOSENSOR

Stick the Biosensor on the chest area, below the collar bone, left of the sternum.

9

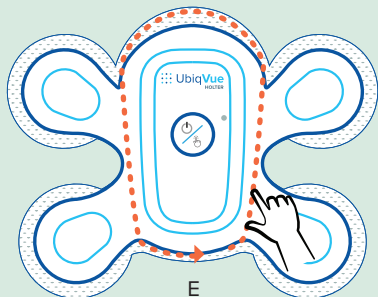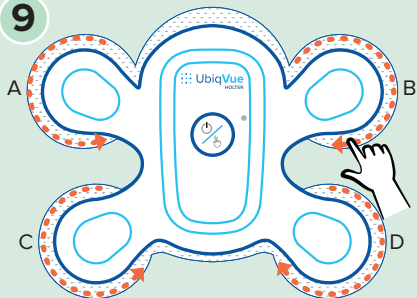

### STICK BIOSENSOR TO SKIN

Press down firmly along the dotted lines with your finger, repeating the steps 10 times.



After Biosensor application, limit your activities for at least one hour so it can fully adhere to your skin.



Keep this guide until the end of the procedure.



## CONNECT THE BIOSENSOR

1



OR





UbiqVue Holter




## DOWNLOAD THE APP

Scan the QR code to download the mobile app on to your smartphone. Or search for UbiqVue Holter in your app library.

2

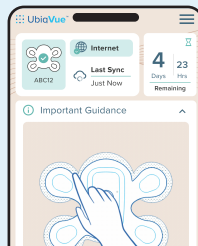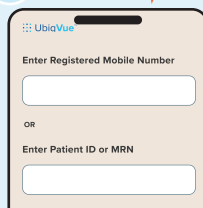

3

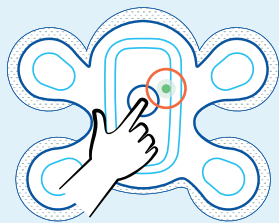


## LOGIN AND TURN THE BIOSENSOR ON


Log in to the UbiqVue Holter App using your mobile number, or patient ID or MRN. Press the Biosensor button— a green light will flash. Follow the app instructions to start the procedure. The procedure is active when there is a green tile in the left corner and your Biosensor ID is displayed.



## USING THE UBIQVUE HOLTER APP

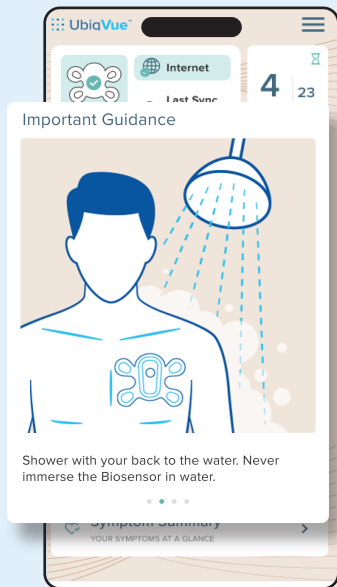

The app will provide important guidance and updates about wearing your Biosensor.

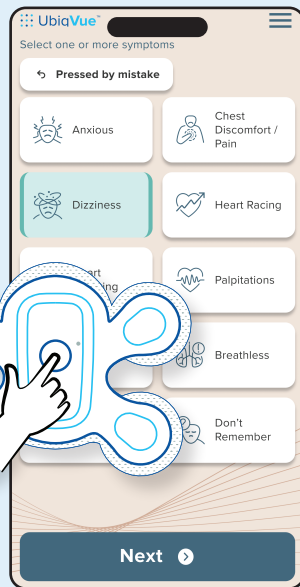


You can record a symptom by pressing the Biosensor button and reporting details in the app.



## BIOSENSOR REMOVAL




**Remove the Biosensor only after the UbiqVue Holter App confirms the procedure is complete.**

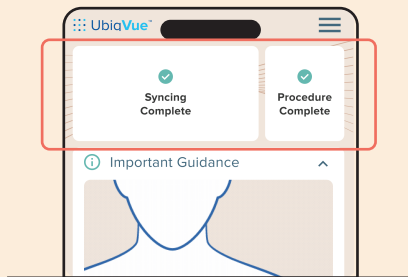

Check the App: Wait until both Syncing Complete and Procedure Complete are displayed.

**1**

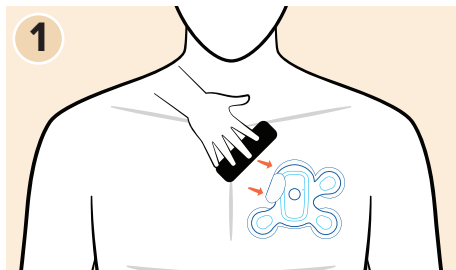


### PEEL OFF THE BIOSENSOR

Slowly peel off the Biosensor starting from a corner. Rub the Remover Wipe on the skin while peeling the Biosensor.



2

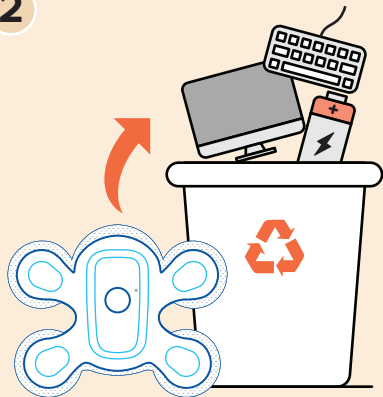

## DISPOSE OF THE BIOSENSOR

Dispose of the Biosensor in accordance with local laws and guidelines for electronic waste.

REF

UC2251-PK



LifeSignals, Inc.  
426 S Hillview Drive  
Milpitas, CA 95035, USA  
[www.lifesignals.com](http://www.lifesignals.com)

EC REP

Renew Health Limited, IDA  
Business & Technology Park  
Garrycastle, Dublin Rd, Athlone  
N37 F786, Ireland

LifeSignals, Inc. © 2025

Doc: 100002117B – February 2025